


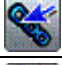


Viewer Toolbar

Icon	Shortcut Keys	Description	Icon	Shortcut Keys	Description
		Close the current study.			Show reports.
		Stop loading the study in progress.		Ctrl + T	Mark study as a teaching study.
		Save the selected (modified) image(s).		Ctrl + K	Open the Sticky Notes window.
	Alt + Space	Mark the selected image as a key image.			Select mode to specify how actions should be applied.
		Open the key image previewer.		Ctrl + I	Put the cursor in neutral mode.
		Open cine mode for creating and playing cine movies. See Cine Toolbar .		W	Apply the windowing mode.
		Export selected images for viewing and printing in the Film Previewer. See Film Previewer Toolbar .		Z	Use the zoom tool.
		Preview the new series.		Ctrl + M	Use the magnifying glass.
		Select the active viewbox layout.		P	Use the pan tool.
		Select the group layout.			Toggle DICOM-defined shutter areas on and off.
		Enlarge the current group.			Flip the selected image horizontally.
		Stack images in the group.			Flip the selected image vertically.
		Manually link images.			Rotate the selected image clockwise.
		Automatically link images.			Rotate the selected image counterclockwise.
		Unlink all groups.		Ctrl + A	Apply automatic windowing to the selected image.
		Preview a report.			Apply inverse windowing to the selected image.
	Ctrl + O	Open the Patient Mini-Archive. See Patient Mini-Archive Toolbar .			Select a graphics tool.